

TESTING YOUR DRUGS FOR FENTANYL

Using BTNX Rapid Response™ Fentanyl Test Strips



WHY SHOULD I TEST MY DRUGS FOR FENTANYL?

- Fentanyl is a powerful opioid. It's often found in pressed pills or powder.
- If you know fentanyl might be in your drugs, it can help you make decisions to lower your risk of overdose.

WHAT SHOULD I KNOW ABOUT FENTANYL TEST STRIPS?

- Fentanyl isn't always mixed evenly into drugs. A single negative test does not mean your entire supply is fentanyl-free.
- Test strips should be used in combination with other harm reduction strategies:
 - Start with a smaller dose or test shot.
 - Take turns using with a friend.
 - Call a friend or the Never Use Alone Hotline (877-696-1996) if you use alone.
 - Carry naloxone and know how to use it.




DOH 150-240, September 2023. To request this document in another format, call 1-800-525-0127. Deaf or hard of hearing customers, please call 711 (Washington Relay) or email doh.information@doh.wa.gov.

HOW TO USE BTNX FENTANYL TEST STRIPS

1. PREPARE YOUR SAMPLE

- To test residue left in a cooker:
 - Add 2mL of clean water (2 full syringes) to your cooker. Mix well.
- To test pills, crystal, and powders:
 - Crush up your drugs and mix them up.
 - Add 10mg of drugs (about the size of a match head or grain of rice) in a cooker.
 - Add 5mL of water (a full cooker) for every 10mg of drugs. Mix well.

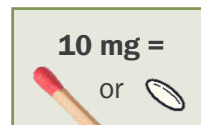

2. DIP THE STRIP FOR 15 SECONDS.

- Hold the dark blue end of the test strip. Don't dip it past the wavy blue lines.

3. WAIT 5 MINUTES FOR RESULTS.

READING RESULTS

One red line means the test is positive. You may have fentanyl in your drugs.

1 RED LINE = POSITIVE for FENTANYL



Two red lines means the test is negative. Look closely for a second line – it can be very faint!

2 RED LINES = NEGATIVE for FENTANYL



! If there is only a red line on the left or no red lines appear, the test is **INVALID**. Redo the test with a new strip.